

# Central Oregon Emergency Mask Makers Original Pleated Mask (Nose Wire + Pocket)

This pattern is excellent for sewing masks for yourself and friends. If you'd like to make masks for St. Charles Health System's 10,000 mask campaign, join the Facebook group, Central Oregon Emergency Mask Makers for all the specifics.



← 7 1/2" (19cm) →

**Central Oregon Emergency Mask Makers**  
4 April 2020 - Page 1/4


**Notes:**

- 100% Cotton, high thread count, tight weave (pillow cases, sheets, quilting fabric)
- Cotton thread please
- Please pre-wash and dry fabric before sewing
- Straps can be cotton bias tape or same fabric
- Light colors are preferred but patterns can be fun
- Please no offensive or divisive symbols/patterns! Keep it fun for all!
- No Elastic Straps (Most elastic has latex or no guarantee of being latex free)
- Avoid pinning if possible, use seam clips or binder clips to hold fabric
- This pattern complies with the St. Charles Medical Center request.


**1 Cut Fabric Body:** 8 3/4" (22cm) Tall x 7 1/2" (19cm) Wide, quantity 2 pieces per mask. \*\*Two different patterns per mask\*\*

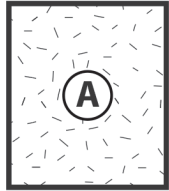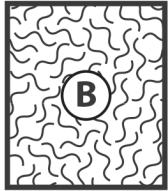**2 Double Fold 1/4" & Top Stitch Bottom Edge:** A & B  
Follow: 2a 2b 2c



This is a full size pattern, check print size here with a nickel or measure 1" (25mm) Square.

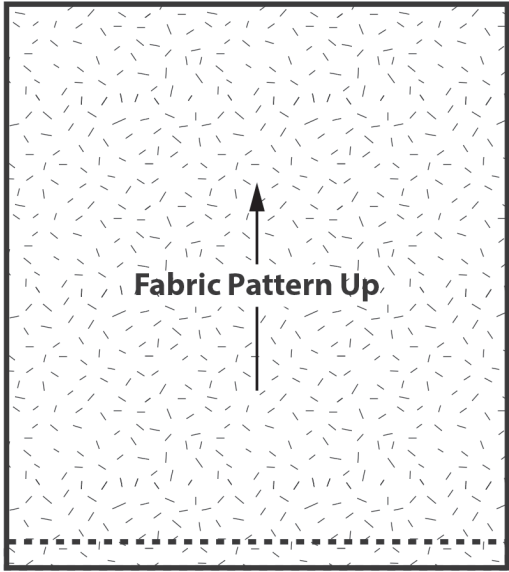


**2a** Fold Under 1/4" (6mm)


**2b** Fold Under 1/4" (6mm)

**2c** Top Stitch Scant 1/4 (5mm) or 1/8" (3mm)

**3 Pleat A & B:** Fold 2 pleats using lines as shown. Shingle down. Iron. Finished size should be 4 1/4" (11cm) tall.  
Follow: 3a 3b Check with: 3c




**3c** 4 1/4" (11cm)

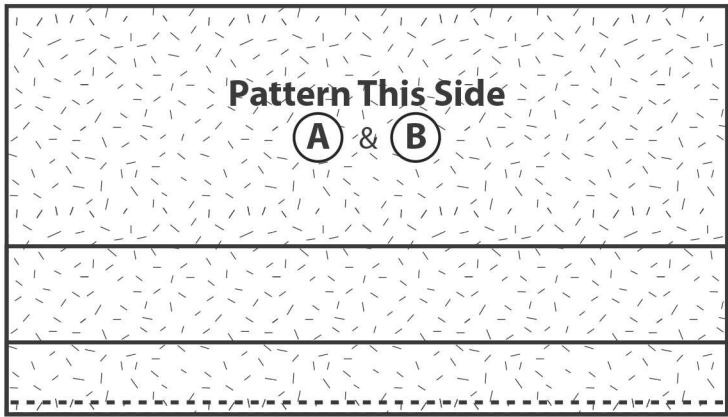

**3a** Fold Over 1" (25mm)

**3a** Fold Over 1" (25mm)

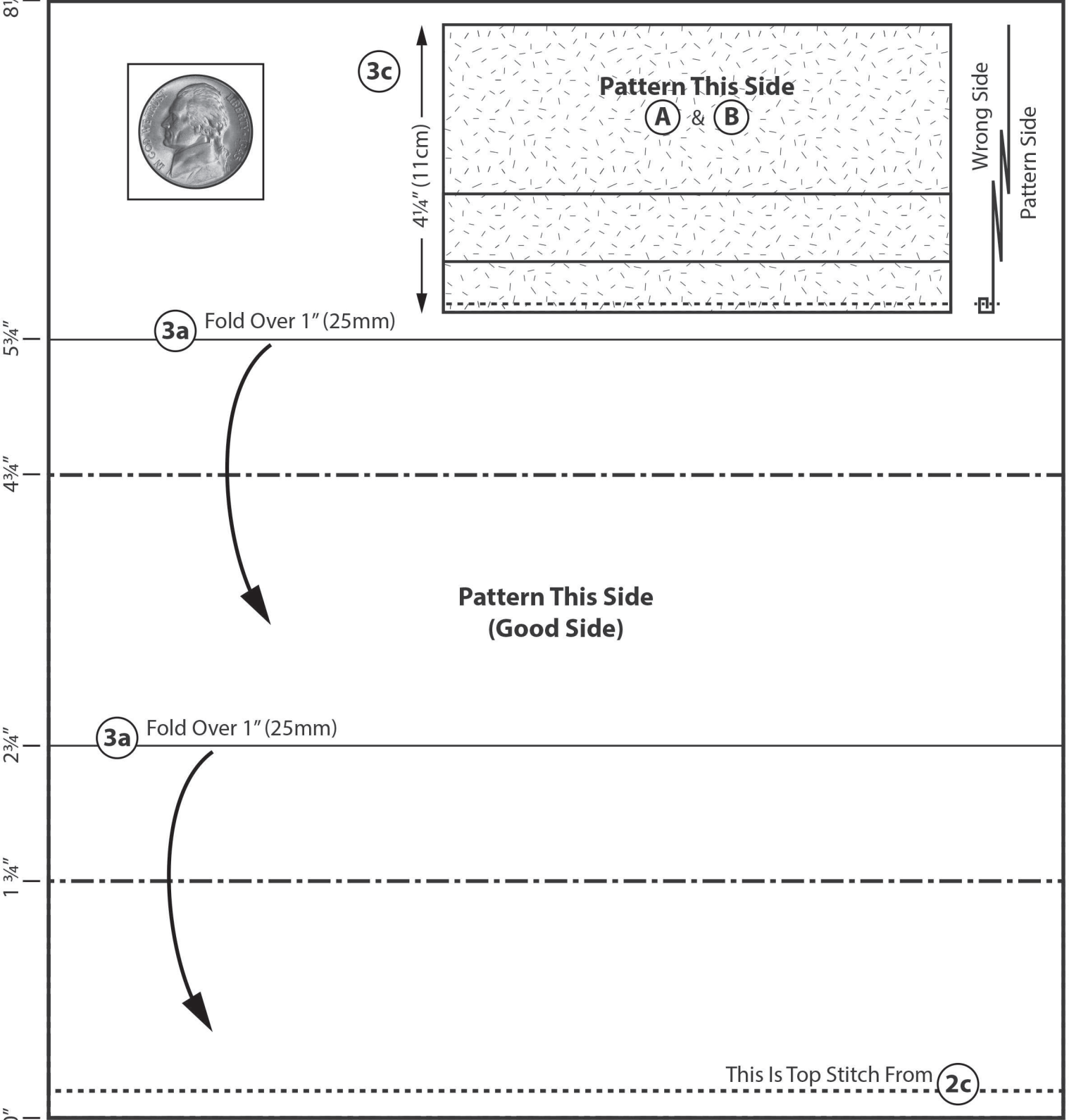

This Is Top Stitch From **2c**